


**eso**

**EHR Cardiac Monitor  
Cloud Tutorial**

# Index

Import a Case .....	<b>3,4</b>
Review a Case .....	<b>5</b>
Un-Import a Case .....	<b>6</b>
Error Messages .....	<b>7</b>

1. From the Vitals Tab, click Monitor Import.

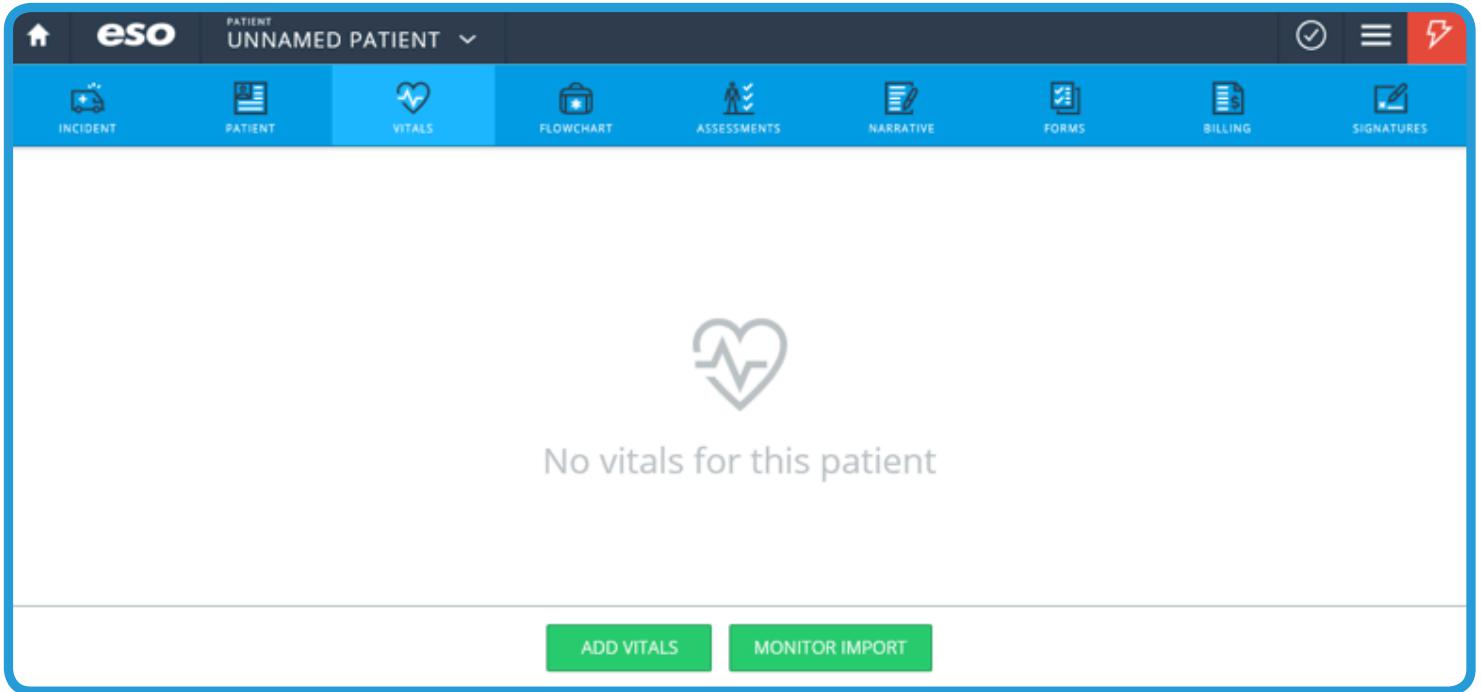

2. From the list of available cases retrieved from the cloud service, select the appropriate case and click Import

**Cardiac Monitor - Select data to import**


PHYSIO CONTROL

Device: All

Start Time	Duration	Device	ID	Patient
03/24/2008 @ 11:03	28m	001 - LP12	85	Quest, Coach
11/29/2007 @ 19:11	52m	006 - LP12	84	Petri,
07/03/2007 @ 9:07	12m	006 - LP12	83	Clicquot, Philippe
03/06/2003 @ 8:03	16m	005 - LP12	86	Zimmer,
08/14/1998 @ 13:08	12h 22m	006 - LP12	79	Murray, SALLY


CANCEL IMPORT

3. The list of available cases can be filtered by selecting a specific device from the **Device** list.

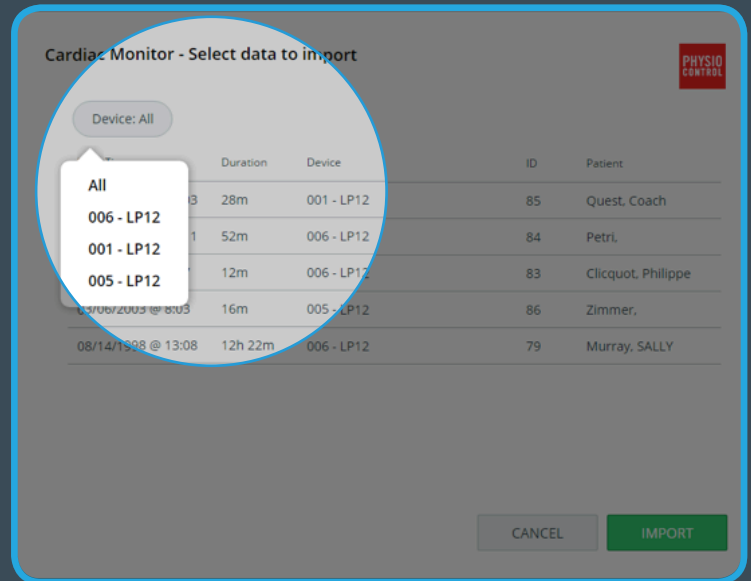

4. When the process is complete, a summary of how many vitals, ECG images, and flowchart interventions were imported into the record. Click **OK** to continue.

**Note:** Imported data should not be used in conjunction with patient diagnosis or provided to a hospital prior to patient arrival.

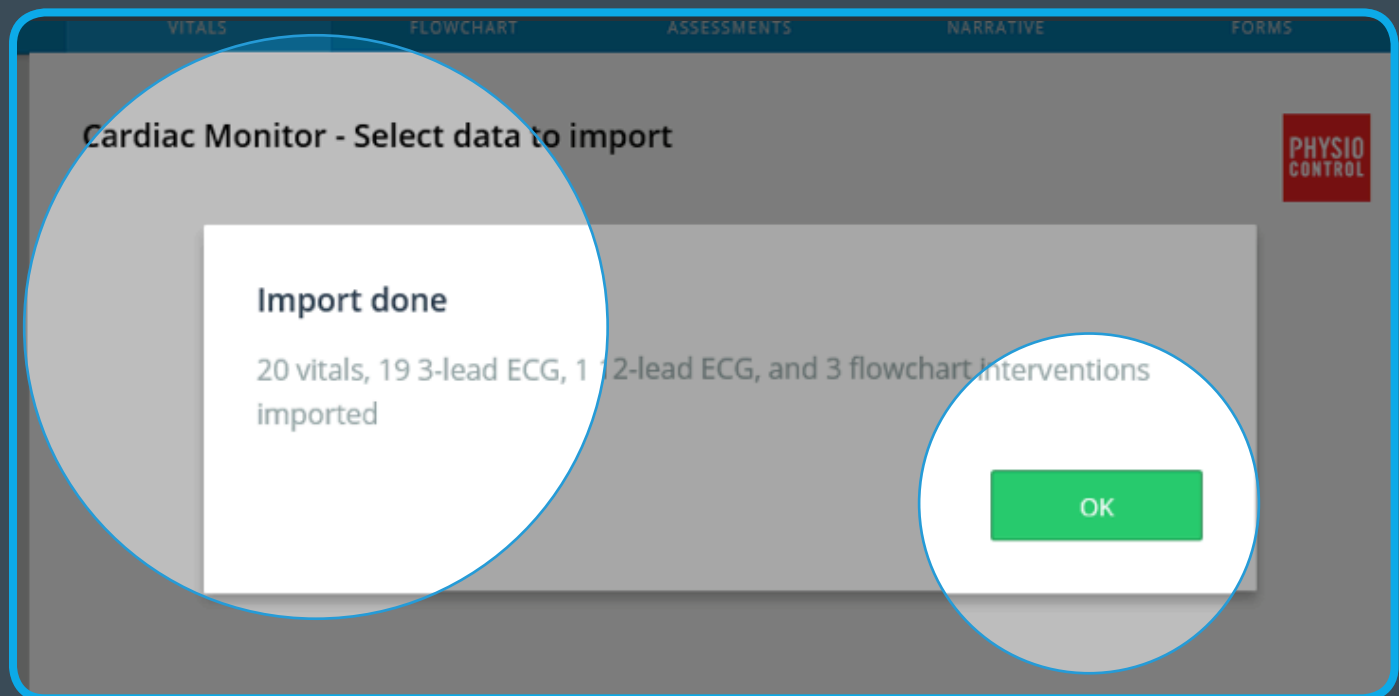

- Each vitals reading is timestamped and broken down into 3 sections of information:

**Note:** The heart icon represents an associated ECG strip. The 3 represents a single lead view, where the 12 represents an image of a 12-lead ECG acquisition.

### General













AVPU, Side, Position, Blood Pressure, MAP, Pulse, Respiration Rate, SpO2, EtCO2, CO, Blood Glucose, Temperature

### Scoring


Revised Trauma, Pediatric Trauma, Glasgow Coma, Glasgow Coma Qualifier, Pain

### ECG

3-Lead details, 12-lead details

<span>General</span> <span>Scoring</span> <span>ECG</span>													
Time	AVPU	Side	Pos	BP	MAP	Pulse	RR	SpO <sub>2</sub>	EtCO <sub>2</sub>	CO	BG	Temp	
 11:54:46	--	--	--		--	65	--		6	--	--		
 11:54:08	--	--	--		--	75	9		8	--	--		
11:53:30	--	--	--		--	76	11		11	--	--		
 11:50:01	--	--	--		--	--	--		--	--	--		
11:48:30	--	--	--		--	132	12		38	--	--		
 11:46:10	--	--	--		--	124	13		38	--	--		
 11:46:05	--	--	--		--	61	13		38	--	--		

- To view the ECG image, select **ECG** from the Row **Menu Icon** on the far right.

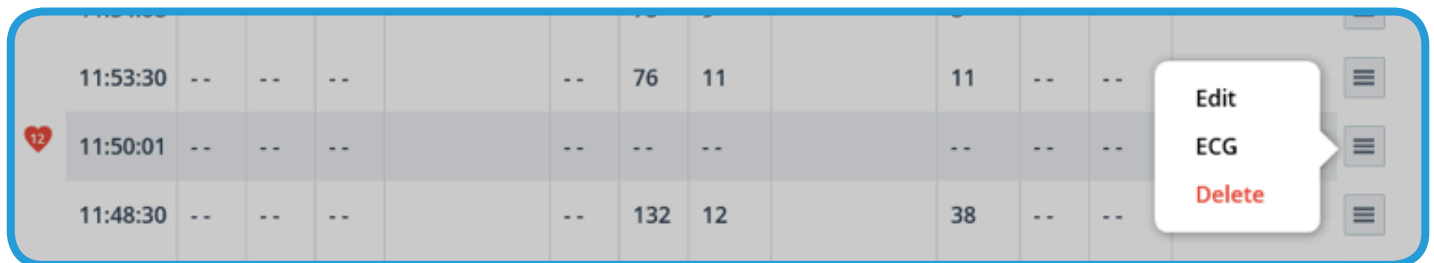

If an ECG file has been wrongly imported, simply select the [Undo Import] option in the bottom right corner.

If more than one file has been imported, **Undo Import** will remove the most recently imported file. You can un-import all files that have been imported by continually selecting **Undo Import** and confirming.

The screenshot shows a software interface with a top navigation bar containing icons for INCIDENT, PATIENT, VITALS, FLOWCHART, ASSESSMENTS, NARRATIVE, FORMS, BILLING, and SIGNATURES. Below the navigation bar are three tabs: General (selected), Scoring, and ECG. A table displays vital signs data with columns for Time, AVPU, Side, Pos, BP, MAP, Pulse, RR, SpO<sub>2</sub>, EtCO<sub>2</sub>, CO, BG, and Temp. The table contains seven rows of data, each with a red heart icon in the Time column. To the right of the table is a vertical menu with three-line icons. Below the table are two green buttons: 'ADD VITALS' and 'MONITOR IMPORT'. At the bottom right, there is a grey button labeled 'Undo Import' with a red arrow pointing to it from the table area.

Time	AVPU	Side	Pos	BP	MAP	Pulse	RR	SpO <sub>2</sub>	EtCO <sub>2</sub>	CO	BG	Temp
11:54:46	--	--	--		--	65	--		6	--	--	
11:54:08	--	--	--		--	75	9		8	--	--	
11:53:30	--	--	--		--	76	11		11	--	--	
11:50:01	--	--	--		--	--	--		--	--	--	
11:48:30	--	--	--		--	132	12		38	--	--	
11:46:10	--	--	--		--	124	13		38	--	--	
11:46:05	--	--	--		--	61	13		38	--	--	

The screenshot shows a confirmation dialog box with a blue border and a blue exclamation mark icon in a circle on the top left. The dialog box has a white background and contains the following text: 'Undo cardiac monitor import?' followed by 'This will remove all vitals and flowchart entries that were associated with this import.' At the bottom of the dialog box are two buttons: a grey 'Cancel' button and a green 'Undo import' button. The dialog box is overlaid on a blurred background of the software interface, showing parts of the vital signs table.

If you see an unable to Connect Error, confirm that you have internet connectivity. If not, wait until you are connected and try again.

If you have connectivity, contact your agency administrator to confirm the web set up is working correctly.



Cardiac Monitor - Select data to import

PHYSIO CONTROL

Device: All

Unable to connect

CANCEL TRY AGAIN